


FORCE LIFTUP ACL/PCL FIXATION DEVICE WITH LOOP

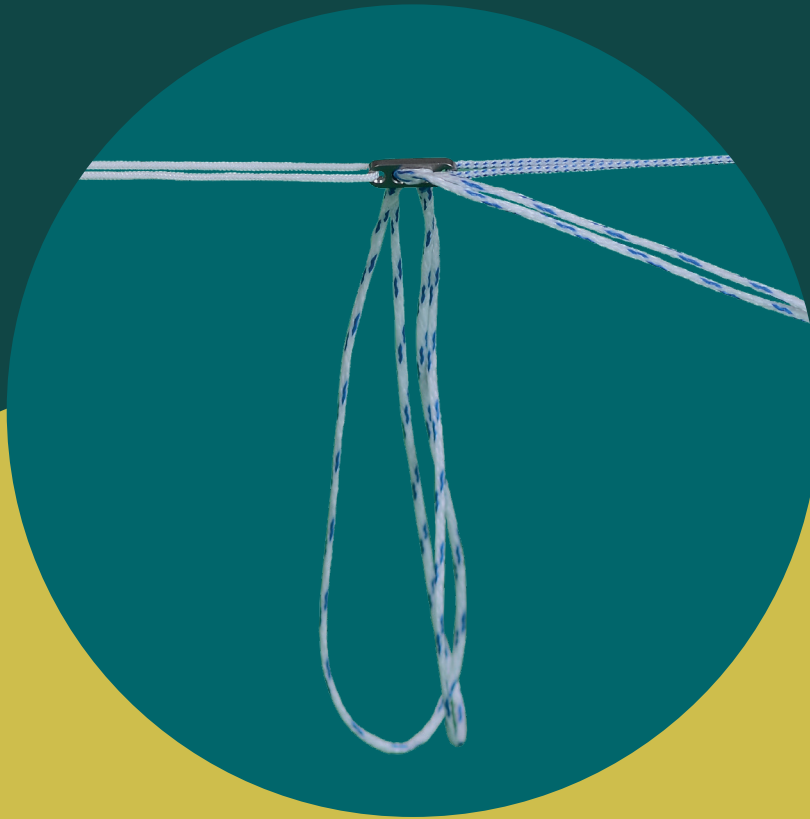

This products are manufacturing by Titanium alloy and UHMWPE fibers for ACL and PCL repair. Button measurements are; lenght 12mm widht 3,90mm and thickness 1.5mm. They can be use with all brand acl instruments with 4.5mm or 5mm cannulated drills. Loop lenghts are adjustable that what the surgeons needs. Loop ropes width is USP #7, and adjusting ropes are USP # 2. This implant's loop lenghts adjusting and locking from femoral side with free rope ends.

"POWERFUL, FAST AND RELIABLE"

www.aledamakina.com



REF CODE

PRODUCT NAME

FORCE PLUS LIFTUP ACL/PCL FIXATION DEVICE WITH LOOP

YN00234

FORCE PLUS LIFTUP ACL/PCL FIXATION DEVICE WITH LOOP

Çamlıca Mah. Anadolu Bulvarı R Blok Apt. No: 20 R/6
Yenimahalle | ANKARA | TURKEY

Phone: +90 312 397 70 41 | **Fax:** +90 312 397 70 42

info@aledamakina.com